

INSTALLATION



2. PLUG IN MODULAR CABLE(#5775) TO ZONE VALVE AND DESK OR WALL SENTINEL ALARM.

3. <u>FUNCTION TEST SYSTEM</u> BY VARYING PRESSURE AND MONITORING ALARM FOR BOTH N₂O OR O₂ AT HI, NORMAL, AND LOW PRESSURE SIGNAL CONDITIONS.

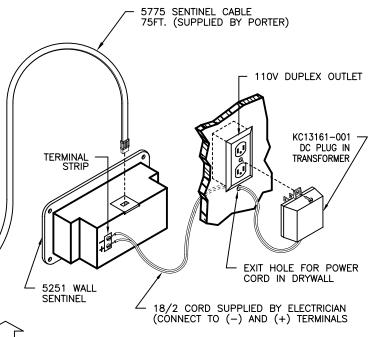

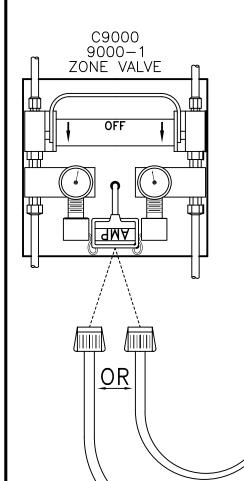

WALL AND DESK SENTINEL WIRING ILLUSTRATION A-2522-000 MODULAR PLUG JUNCTION BOX ADDITION BOX

5201 DESK
SENTINEL

C-637-000

DC PLUG-IN TRANSFORMER
WITH 6FT. POLARIZED LINE
CORD TO 110V AC SUPPLY.

5775 SENTINEL CABLE, 75FT
(SUPPLIED BY PORTER)

PROPRIETARY NOTICE
THIS DOCUMENT AND THE DATA DISCLOSED HEREIN IS
NOT TO BE REPRODUCED, USED OR DISCLOSED WHOLE
OR IN PART TO ANYONE WITHOUT THE PERMISSION OF
PORTER INSTRUMENT DIMISION, PARKER HANNIFIN CORPORATION.

DATE MAR. 13, 1986 DR. ALE

REV EC NO. DESCRIPTION

APP. JDD

NAME | DATE

APPVD

REV

REVISIONS

SCALE NONE

Parker Hannifin Corporation Porter Instrument Division 245 Township Line Rd, PO Box 907 Hatfield, PA 19440-0907, USA

CH. JDD

WIRING DETAILS FOR ZONE VALVE, MODEL 9000-1 AND C9000 **FM-361** D SH 1 OF 1

FILENAME: 0035110D

NUMBER

PORTER